

**ICX-M2**  
Intensive Care Bed



3/4 Length  
Siderails

Sturdy Linkages

Full Set of  
Clinical positions




*#polymerheadandfootpanels #standardconfiguration*

# ICX-M2

## Intensive Care Bed

### Clinical Positions

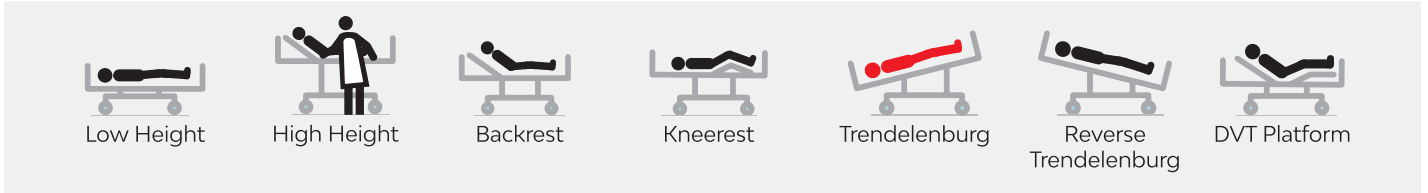

### Features

- Backrest, kneerest, height, Trendelenburg and reverse Trendelenburg adjustable on crank mechanism.
- Four section perforated mattress platform.
- The adjustable calf rest adapts to obtain a good DVT platform.
- The caregiver can obtain full set of clinical positions in the ICX bed.
- Anti-microbial powder coating.

### Specifications

Model	ICX-M2
External dimension (L X W) mm	2190 X 1020
Mattress platform dimension (L X W) mm	1965 X 860
Height adjustment (H) mm (without mattress)	525 to 790
Backrest angle	0° to 70°
Kneerest angle	0° to 35°
Trendelenburg angle	0° to 10°
Reverse Trendelenburg angle	0° to 7°
SWL - Safe working load	135 kgs
Patient weight bearing capacity	110 kgs

Engineering variation  $\pm 2^\circ$ ,  $\pm 10$ mm

### Standard Accessories

Polymer moulded head and foot panels (PP2) - easily removable
3/4 Length siderails
125 mm $\varnothing$ high grade synthetic body casters with brake on two diagonally opposite casters (C3)
Removable chromium plated forged handle
Buffers on all four corners
IV pole provision
SS 304 heavy duty IV pole with two hooks
In-built urine bag holder

### Optional Accessories

Four section 40 density foam mattress of 100 mm thickness covered with rexine (M1)
3/4 Length Vertica siderails
MS oxygen cylinder cage
SS 304 telescopic IV pole with two hooks
Removable polymer moulded handle

Subject to technical changes without notice. Accessories in the picture may not be part of the standard configuration. The product colours in the pictures may differ slightly from the actual product.

#### MUMBAI OFFICE

AMSON BIOMED Pvt. Ltd.,

1st Floor, Lal Bungalow, Upashray Lane, Ghatkopar East, Mumbai 400077

Phone +91 9820032799 | info@amsonbiomed.com | www.amsonbiomed.com

