

## >> Transcure-P

Trauma & Emergency Care Recovery  
Stretcher Trolley



Refreshing  
New Look

Sturdy  
Hydraulic Pump

Ergonomically  
Designed Siderails




*#safeexpress*

# Transcure-P

## Trauma & Emergency Care Recovery Stretcher Trolley

### Clinical Positions

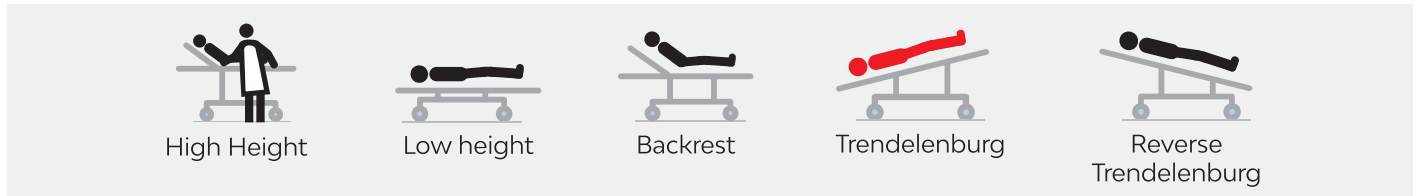

### Features

- Height adjustment by foot operated hydraulic pump.
- Trendelenburg and reverse Trendelenburg operated on gas spring.
- Backrest adjustable on ratchet.
- X-ray permeable removable two section stretcher top made of HPL (high pressure laminate).
- Ergonomically designed handles for better manoeuvrability.
- Ergonomically designed SS 304 swing away type siderails.
- Anti-microbial powder coating.

### Specifications

Model	Transcure-P
External dimension (L X W) mm	1985 X 785
Mattress platform dimension (L X W) mm	1830 X 630
Height adjustment (H) mm (without mattress)	690 to 940
Backrest angle	0° - 65°
Trendelenburg angle	0° - 17°
Reverse Trendelenburg angle	0° - 6°
SWL - Safe working load	165 kgs
Patient weight bearing capacity	135 kgs

Engineering variation  $\pm 2^\circ$ ,  $\pm 10$ mm

### Standard Accessories

125 mm $\varnothing$ high grade synthetic body casters with brake on two diagonally opposite casters (C3)
Mattress of 50 mm foam with rexine
Sliding X-ray cassette holder
Safety belts
SS 304 telescopic IV pole with two hooks
MS oxygen cylinder cage
Buffers on all four corners

### Optional Accessories

Storage facility at both ends
SS 304 heavy duty IV pole with two hooks

Subject to technical changes without notice. Accessories in the picture may not be part of the standard configuration. The product colours in the pictures may differ slightly from the actual product.

#### MUMBAI OFFICE

AMSON BIOMED Pvt. Ltd.,

1st Floor, Lal Bungalow, Upashray Lane, Ghatkopar East, Mumbai 400077

Phone +91 9820032799 | info@amsonbiomed.com | www.amsonbiomed.com

